



SKB041-3

PATIENT TRANSPORTATION TROLLEY

Registration code: SKM-G




SAIKANG / Medical

Healthier
your life

SKB041-3

FEATURES AND HIGHLIGHTS

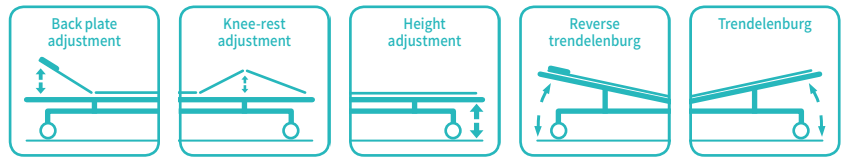

Technical parameters

X-ray Backrest

▪ External size (LxW):	2120x840 mm
▪ Height from floor to bed platform:	550-860 mm
▪ Height from floor to lowered guardrail:	590-900 mm

X-ray Platform

▪ External size (LxW):	2120x840 mm
▪ Height from floor to bed platform:	600-910 mm
▪ Height from floor to lowered guardrail:	640-950 mm

▪ Back plate adjustment:	0-85° (±5°)
▪ Knee-rest adjustment:	0-40° (±5°)
▪ Reverse Trendelenburg:	0-15°
▪ Trendelenburg:	0-15°



Technical configuration

<input checked="" type="checkbox"/> Usa hydraulic pump	2 pcs
<input checked="" type="checkbox"/> Crank system	1 set
<input checked="" type="checkbox"/> Back rest air spring	2 pcs
<input checked="" type="checkbox"/> 8" German castors,central locking	4 pcs
<input checked="" type="checkbox"/> Emergency push handle	1 set
<input checked="" type="checkbox"/> Directional wheel	1 pc
<input checked="" type="checkbox"/> Folding side rail	1 set
<input checked="" type="checkbox"/> Φ98mm wall bumper wheel	4 pcs
<input checked="" type="checkbox"/> Iv pole previsions	4 pcs
<input checked="" type="checkbox"/> Mattress	1 pc
<input checked="" type="checkbox"/> IV pole	1 set
<input checked="" type="checkbox"/> Central locking break pedal	2 sets
<input checked="" type="checkbox"/> Directional wheel release pedal	2 pcs
<input checked="" type="checkbox"/> X-ray translucent bed platform	1 set



Metal material

Baosteel from Fortune Global 500
Panasonic Robotic ensure 360° full
smooth wedding.

**Plastic**

Environmental plastic can be
sterilized in 100°C, tensile upto
30MPa.

**Steels | Painting**

11 Process epoxy painting, ASTM
testing anti-bacterial, paint thickness
0.12 mm, brightness 60°, paint can
resist 50 kg impact.

**Brake**

Break spring keep qualified after
10,000times uninteruped break.

**Hydraulic pumps**


Height adjustment, trendelenburg
and revised tendelenburg through
two USA hydraulic pumps.

**Center wheel**

Can be engaged to facilitate direction
control, enabling users to navigate
the stretcher easily in hectic and
dynamic environments.

**Rising**

Hydraulic treadle for rise.

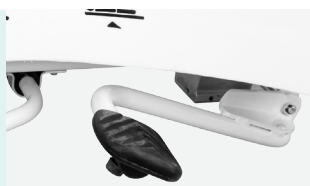**Falling**

① Leg section lows down to make
head higher.
② Back section lows down to make
leg section higher.
③ Step on the middle of the pedal,
make the bed base fall down steadily.

**Castors**

Four double side casters of 8 inch
diameter, simultaneously lockable by
the step pedal at both sides.
TPR tire no worn out after running
30KM, save anti-winding hard shell,
united forming without bolts.

**Mattress**

Waterproof PU mattress, convenient
for transfer patient.

**Side rail**

Foldable aluminium alloy guardrail

**Bumper wheel**


Bumper wheel detachable, max.
protect the bed and patient during
transportation

**Drainage Hooks**

Mobile plastic drainage hooks.


**IV pole**

Stainless Steel Column with plastic
hook.

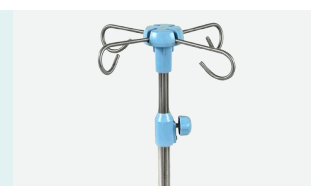



CONSTANT INNOVATION

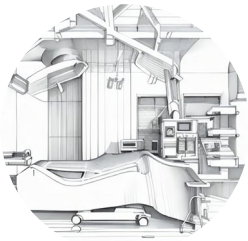


Design



Research



Manufacturing

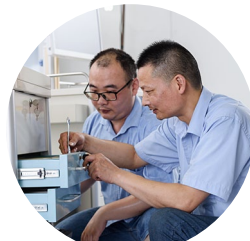

Maintenance



Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc
W3400 Corrientes, Argentina
TEL: +55 (11) 2758-0438
Rua Universal, Hollywood,
São Bernardo do Campo,
Brazil
TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd
Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

Saikang Medical Technology Ltd
Sala Thammasop Subdistrict
Thawi
Watthana District
Bangkok Thailand
TEL: +66 642741046

UK Branch

Uk Healthward International Limited
DEPT 111,
196 High Road Wood Green,
London, N22 8HH
United Kingdom
TEL: +31 (0) 202111106

Russia Branch

LLC 'SAIKANG'
454004, Russian Federation,
Chelyabinsk,
Universitetskaya
Naberezhnaya, Russia
TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC
Musaffah
P.O Box 53200, Abu Dhabi,
UAE
Tel: +971 585710652

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue