



# SKE941

RECLINER CHAIRS FOR  
ELDERLY



**SAIKANG** / Medical

Healthier  
your life

# SKE941

## FEATURES AND HIGHLIGHTS

### Technical parameters

▪ External size	1100×770×1280mm
▪ Seat size	540×460×560mm
▪ Tray size	480×255mm
▪ Back-rest adjustment	37°-61°
▪ Leg-rest adjustment	13°-36°
▪ Safe working load	113kg



#### Design

Manual recliner with four sections (back, hip, leg, foot), ergonomic design, ideal for dialysis, blood collection, respiratory and cardiac care, and adapts to patient needs.



#### Mechanical Function

Pneumatic adjustment of the backrest, three modes: sitting, semi-lying, lying; footrest is completely hidden when the chair is in the sitting mode, to prevent bumps and falls.



#### Surface material

Using shaped foam sponge, the thickness is up to 130mm, 1.2mm PU leather cover, fire retardant, handmade finishing; sturdy steel frame for easy cleaning and maintenance.  
Option color: red/ green/ blue.




#### Armrest

Ergonomically widened armrests, with a width of 115mm, providing a safe and comfortable environment for the patient; PU leather cover, which is comfortable and soft.



#### Caster

5-inch caster with brake, easy to drive, the loading capacity up to 113kg.

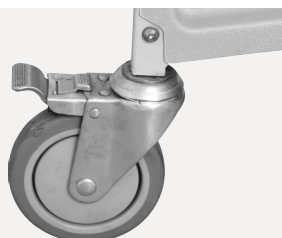

#### Tray

The tray adopts the pull-out installation, and has anti-slip raised edges to prevent teacups and tableware from slipping out; trays on both sides, increased capacity of items to place, easier to pick and place items.



#### Armrest

With armrests on both sides to facilitate the natural stretch of the arms.





## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### South America Branch

**Saikang Medical Technology, Inc**  
W3400 Corrientes, Argentina  
TEL: +55 (11) 2758-0438  
Rua Universal, Hollywood,  
São Bernardo do Campo,  
Brazil  
TEL: +55 (11) 2758-0438

### East Africa Branch

**Saikang Medical Technology Ltd**  
Mezzanine Floor, Taj Tower,  
upper hill, nairobi, Kenya  
TEL: +254714084806

### South East Asia Branch

**Saikang Medical Technology Ltd**  
Sala Thammassop Subdistrict  
Thawi  
Watthana District  
Bangkok Thailand  
TEL: +66 642741046

### UK Branch

**Uk Healthward International Limited**  
DEPT 111,  
196 High Road Wood Green,  
London, N22 8HH  
United Kingdom  
TEL: +31 (0) 202111106

### Russia Branch

**LLC 'SAIKANG'**  
454004, Russian Federation,  
Chelyabinsk,  
Universitetskaya Naberezhnaya,  
Russia  
TEL: +7 (919)300-00-06

### Mid East Branch

**Saikang Medical Supplies (Mid East) LLC**  
Musaffah  
P.O Box 53200, Abu Dhabi,  
UAE  
Tel: +971 585710652

# SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue