



# SKS036

BEDSIDE TABLE



**SAIKANG** / Medical

Healthier  
your life

# SKS036

## FEATURES AND HIGHLIGHTS

### Technical parameters

- Size (LxWxH) 440×440×790 mm

### Technical configuration

- Castor 4pcs



#### Material

The cabinet is made of first-grade cold-rolled steel plate, the countertop is made of stainless steel 304, the thickness is 0.7mm, the drawer panel and door panel are 0.7mm, and the rest is 0.6mm.




#### Track&Key

It adopts full-scale steel ball mute slide rails, and there is no abnormality in 10,000 push-pull tests. Each drawer bears 35KG and is equipped with a folding key.




#### Table top

Stainless steel material, clean and bright appearance, easy to clean, and strong corrosion resistance.

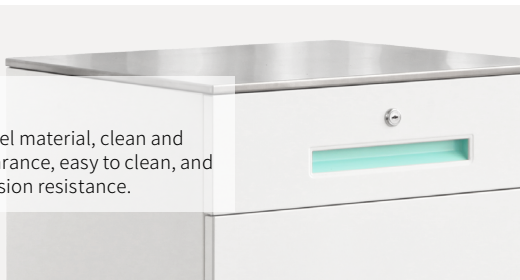


#### Hinges&Magnetic suction

The 304 hydraulic buffer hinge is durable, and the cabinet door is equipped with two magnetic suction on the upper and lower sides, which can realize the gentle closing of the cabinet door and protect the cabinet door.

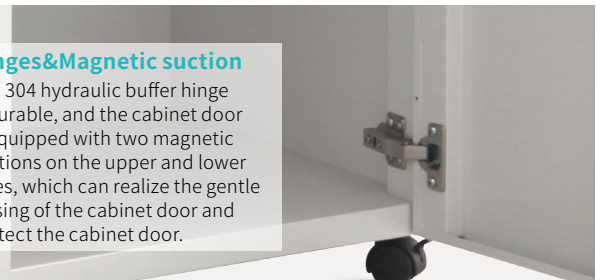

#### Buckle hand

The recessed handle is ergonomic and shows vitality.

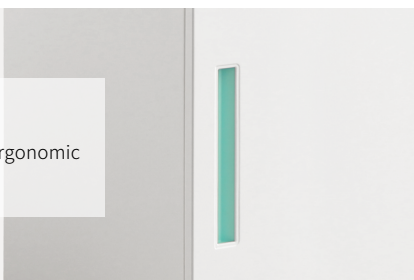

#### Castors

Universal wheel design, the front 2 wheels have their own brakes, it is not easy to move after locking, can travel 5 kilometers, and bear a load of 100KG.





## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### South America Branch

**Saikang Medical Technology, Inc**  
W3400 Corrientes, Argentina  
TEL: +55 (11) 2758-0438  
Rua Universal, Hollywood,  
São Bernardo do Campo,  
Brazil  
TEL: +55 (11) 2758-0438

### East Africa Branch

**Saikang Medical Technology Ltd**  
Mezzanine Floor, Taj Tower,  
upper hill, nairobi, Kenya  
TEL: +254714084806

### South East Asia Branch

**Saikang Medical Technology Ltd**  
Sala Thammasop Subdistrict  
Thawi  
Watthana District  
Bangkok Thailand  
TEL: +66 642741046

### UK Branch

**Uk Healthward International Limited**  
DEPT 111,  
196 High Road Wood Green,  
London, N22 8HH  
United Kingdom  
TEL: +31 (0) 202111106

### Russia Branch

**LLC 'SAIKANG'**  
454004, Russian Federation,  
Chelyabinsk,  
Universitetskaya Naberezhnaya,  
Russia  
TEL: +7 (919)300-00-06

### Mid East Branch

**Saikang Medical Supplies (Mid East) LLC**  
Musaffah  
P.O Box 53200, Abu Dhabi,  
UAE  
Tel: +971 585710652

# SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎️ +86 18021231901 · 🕒 24/7

✉️ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue