



X36

EXAMINATION TABLE

Registration code: SKD-A




SAIKANG / Medical

Healthier
your life

X36

FEATURES AND HIGHLIGHTS

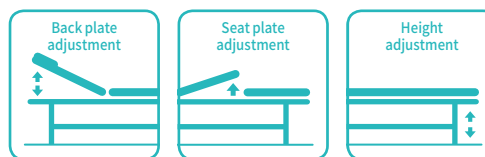

Technical parameters

▪ External Size (LxWxH)	1890x680x470-910 mm
▪ Back Plate Upward Fold	≥80°
▪ Back Plate Downward Fold	≥40°
▪ Leg Plate Upward Fold	≥48°
▪ Leg Plate Downward Fold	≥45°
▪ Bearing Capacity	≥150kg
▪ Headrest Hole	150 mm
▪ Input Voltage	110-230v
▪ Output Power	24W



Technical configuration

<input checked="" type="checkbox"/> Motor	1pc
<input checked="" type="checkbox"/> Pillow	1pc
<input checked="" type="checkbox"/> Gas Spring	2pcs
<input checked="" type="checkbox"/> Brake Pedal	4pcs
<input checked="" type="checkbox"/> 3" Double Wheel	4pcs
<input checked="" type="checkbox"/> Manual Operator	1pc




Optional

<input type="checkbox"/> Bed Surface Foot Lift	1pc
<input type="checkbox"/> Roll Up Tissue Holder	1pc



Bed Frame

The bed frame surface is painted with 1.5mm powder coated steel, and welded by Panasonic Robotic seamlessly. Strong corrosion resistance, environmentally-friendly, bright surface, maintain an overall beautiful appearance. With this structure, the whole bed maintains more stable and makes patients feel more safe.

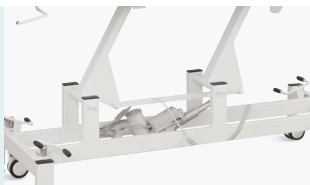

Mattress

The bed platform is integrally made of the high-end composite leather and hand-made cotton mattress, of which the curve design is ergonomic, elastic and soft, breathable, environmentally-friendly, easy to be cleaned, anti-bacterial, impact-resistant.



Pneumatic Pole Support

Two air springs are equipped on the back and legs to provide users with a quiet and safe environment during diagnosis and avoid interference.



Control System

CE certified linear actuator motor, one-button free lift, low noise, and long service life. Manual (foot control) electric overall lifting, manual back-rest lifting, leg-rest lifting adjustment, height adjustment.



Hand Controller

The bed is equipped with Hand-held controller, so that the optimal height is easy to be adjusted by pressing the single button gently. The hook on the back is convenient to take down for patients and nurses.



Bed foot equipment

Equipped with 4pcs 3-inch double-sided wheels, 2 functions side brake with 4 pedals, which are convenient to move and fix. The design of 360° rotating pulley made of polyurethane PU is strong, durable, non-slip, strong in pressure resistance, and good in stability.



Assorted Pillow

The neoteric design of the pillow is narrow on the top and wide on the bottom. Its well-shaped and elegant design matches with modern aesthetic ideas and can blend into the environment harmoniously. Pillow can be disassembled according to demand.



Round Covered Edge

The bed platform was covered in 360°. The soft-touch of surface would be safer for users.




Design

Bed surface two fold three ergonomic design, adjust the height according to the patient's height, personalized design help diagnosis. Dual use of bed and chair, high bearing lasting stability, no fear of heavy pressure impact, to ensure safe use.





CONSTANT INNOVATION

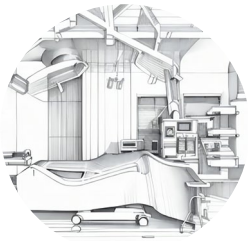


Design



Research



Manufacturing

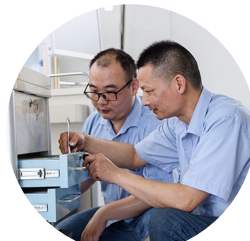

Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue